

AN IMPROVED TECHNIQUE OF COMPUTER- ASSISTED IMAGE ANALYSIS TO QUANTIFY MOLECULAR STRESS IN CETACEANS

¹Anne C. Allen, ^{1,2}Nicholas M. Kellar, ¹Šárka O. Southern and ¹Andrew E. Dizon

¹NMFS, Southwest Fisheries Science Center, La Jolla, California 92038, USA

²Scripps Institution of Oceanography, University of California San Diego, La Jolla, California 92093-02082, USA

Contact email: anne.allen@noaa.gov

Web site: <http://swfsc.nmfs.noaa.gov/>

ABSTRACT

As part of an effort to develop a test for chronic stress in cetaceans, we needed to quantify the expression levels of stress-response proteins (SRP) in immunohistochemically stained skin. We have developed a custom method of computer-assisted image analysis (CIA) using the commercially available software Image-Pro Plus. The immunoreactive signal, corresponding to the SRP expression level, was calculated based on the product of the average staining intensity and the percentage of positively stained epidermal cells. This new method allows for exact measurement of signal directly in the area-of-interest. We demonstrate the application of CIA for quantification of immunohistochemical staining through expression analysis of metallothionein in 22 specimens of beluga, bottlenose dolphins, and spotted dolphins. Image analysis results were compared with the more standard, qualitative method of human-scored visual assessment using light microscopy. Both visual and CIA methods can accurately separate known healthy and diseased cetaceans. However, CIA is able to distinguish staining differences on a much finer scale than a subjective observer can do visually. We conclude that our method of CIA is both accurate and repeatable and provides an efficient means for rapidly quantifying the expression levels of stress-response proteins.

METHODS

- Metallothionein was visualized in thin skin sections from 22 specimens using immunohistochemical (IHC) staining.
- Using identical microscope and camera settings, three digital images per sample were taken to accurately reflect the overall staining. In a separate session, three additional images per sample were taken to assess reproducibility.
- All images were analyzed using the commercially available software program Image-Pro Plus version 4.1 for Windows 95/98/NT (Media Cybernetics, Silver Spring, MD).
- Hardware system used: Gateway E-4200 computer and VX900 monitor, DVC high resolution digital camera 1310C, Olympus System light microscope model BX50 (Scientific Instrument Company, Sunnyvale, CA).

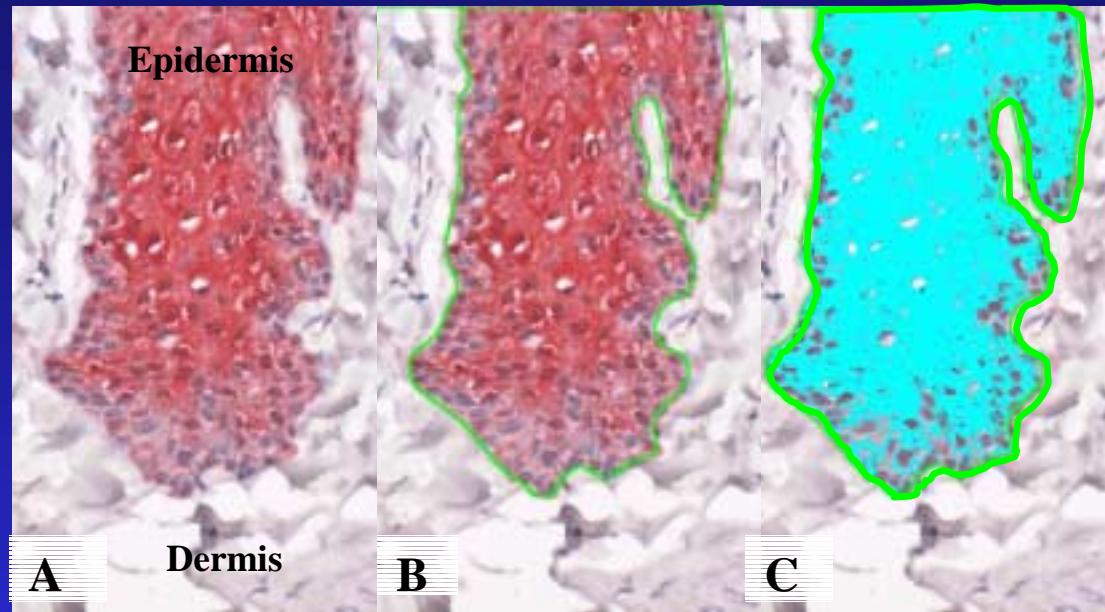


Figure 1. Quantification of the expression of the SRP metallothionein in beluga skin using computer-assisted image analysis.

Metallothionein (MT) was detected in a 5-micron section of beluga skin using immunohistochemical staining with rabbit antiserum to horse MT. As a result of the IHC staining, MT was visualized as a red pigment. MT expression was localized to the immature keratinocytes in the epidermal ridges (panel A). Using the free-form drawing tool in Image-Pro Plus, we selected the epidermis as the area-of-interest. All calculations were performed only within this selected area (panel B). A color file was created that exactly selected the hue, saturation, and intensity that indicated the protein expression. This color file defined the range of the signal and was applied to all samples. When superimposed on the image, the color file was visualized as an artificial blue color that showed the localization of the MT signal (panel C). (Original magnifications x400).

QUANTIFICATION OF THE STRESS-RESPONSE PROTEIN EXPRESSION LEVEL SCORE (ELS)

**ELS = Mean Optical Density of Red x Percent Area
Positively Stained x 100**

- For Figure 1, $ELS = 0.230 \times 0.725 \times 100 = 16.675$
- Observed Range of ELS: 0-35

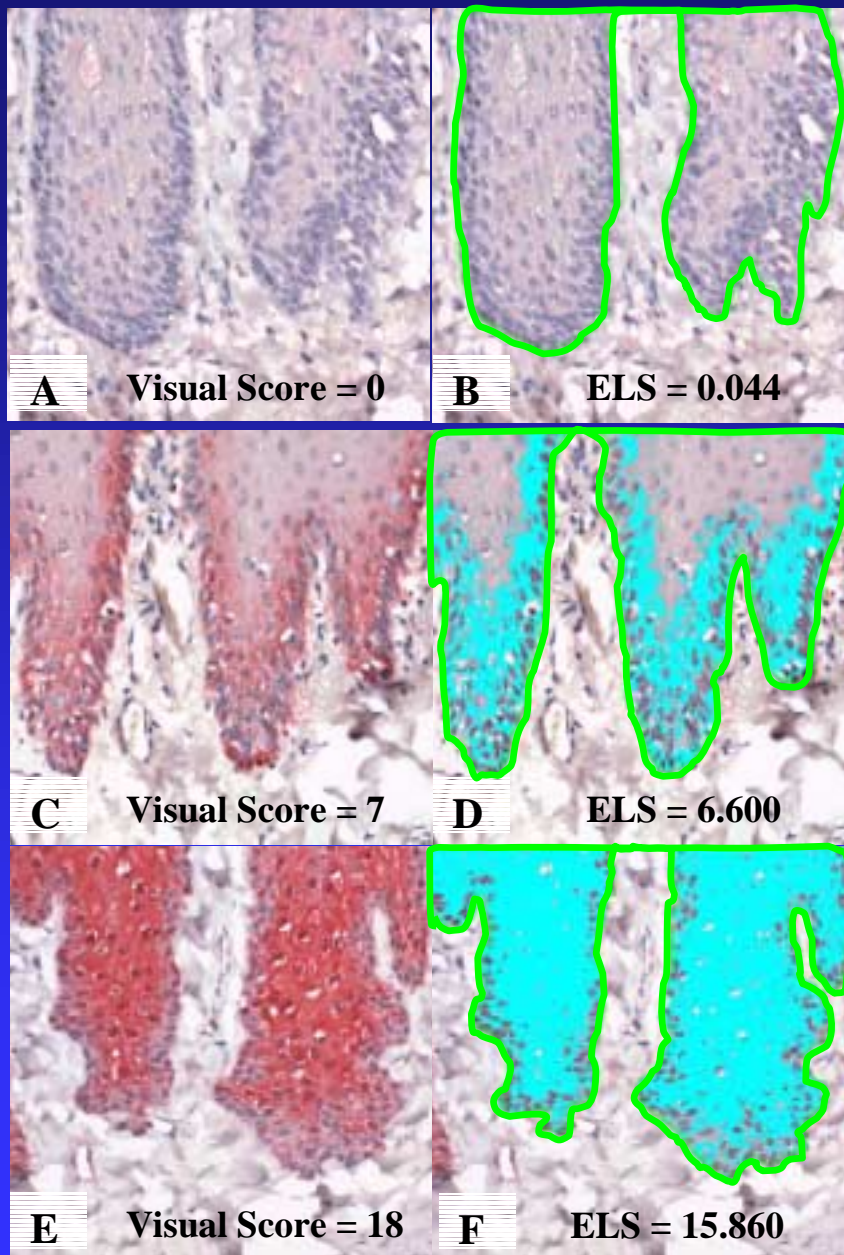


Figure 2. Application of CIA method to cetacean skin sections, comparing with visual analysis.

MT was visualized in 22 different skin specimens, 3 of which are shown here (panels **A** and **B**, **C** and **D**, **E** and **F**). MT levels were evaluated using two methods: subjective visual analysis (panels **A**, **C**, and **E**) and computerized image analysis (panels **B**, **D**, and **F**). Visual analysis involved three repeat evaluations of the specimens using light microscopy at 40x, 100x, and 400x magnifications. CIA was performed on three images per specimen at 400x magnification, and the average ELS score was obtained. Only one image per specimen is shown here; however, the expression level score shown is the average ELS (panel **B**, **D**, and **F**). The ELS found using the computer-assisted method was concordant with the visual assessment of staining (Kendall Coefficient of Concordance, $W = 0.975$, $0.001 < p < 0.01$). (Original magnifications x400).

REPRODUCIBILITY

For all 22 specimens, 3 additional images were taken in a separate session. CIA was performed in the same manner on these images, and again the average ELS value was calculated. The two CIA sessions were shown to give nearly identical ELS values for every specimen (Kendall Coefficient of Concordance, $W=0.993$, $0.001 < p < 0.01$).

See next slide for Figure 3.

Figure 3. Graph of MT expression level scores for both CIA and visual analyses of 11 reference specimens and 11 field specimens.

All 22 specimens were analyzed via subjective visual analysis and computer image analysis. Visual analysis involved three repeat evaluations of the specimens using light microscopy at 40x, 100x, and 400x magnifications. CIA was performed on three images per specimen at 400x magnification, and the average ELS value of all three images was obtained. The eleven field specimens had a range of ELS values from 0-18. Four of the field specimens fell within the range of the known normals, four were within the range of the clinically diseased samples, and three had intermediate ELS values.

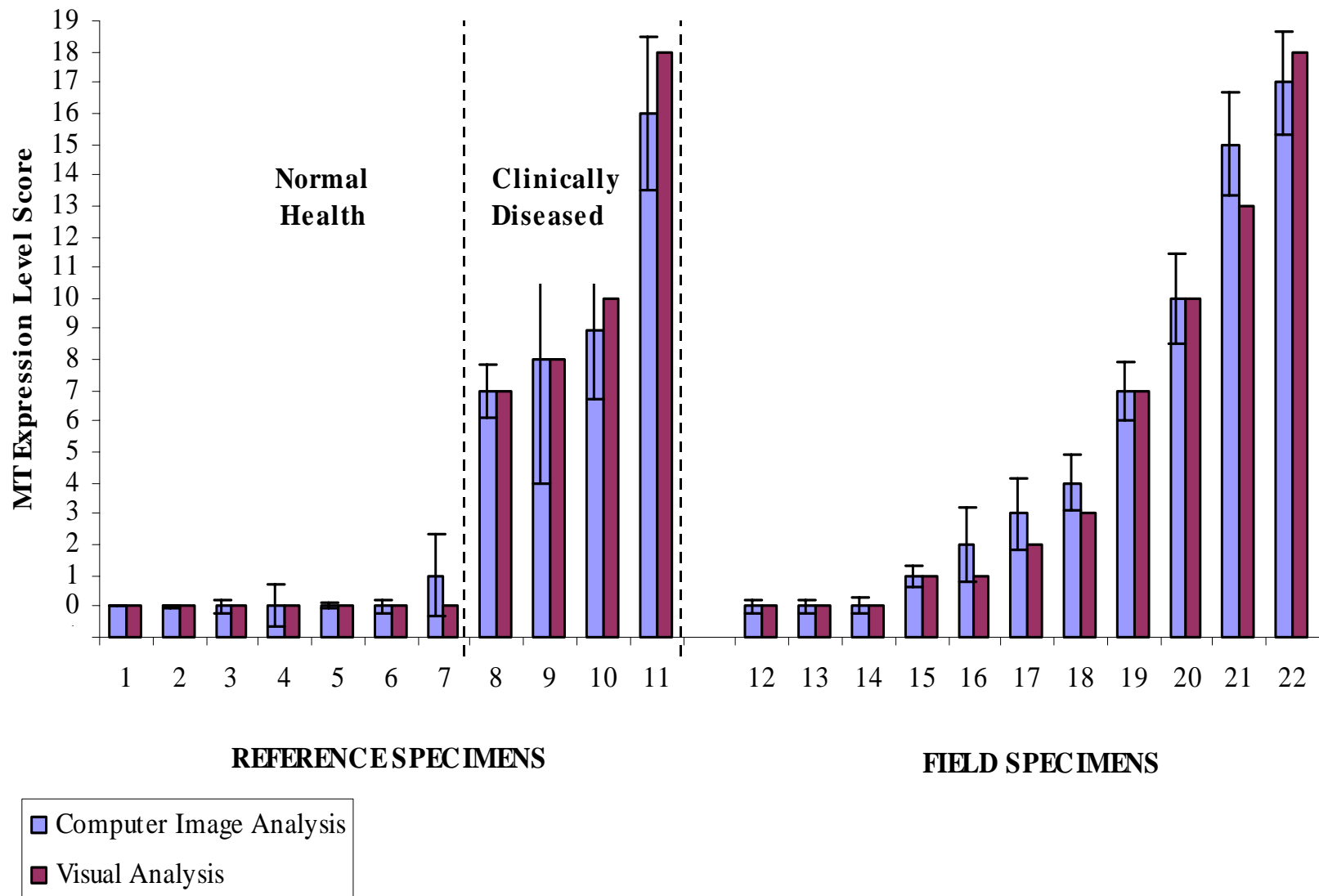


Figure 3. (Description on previous slide)

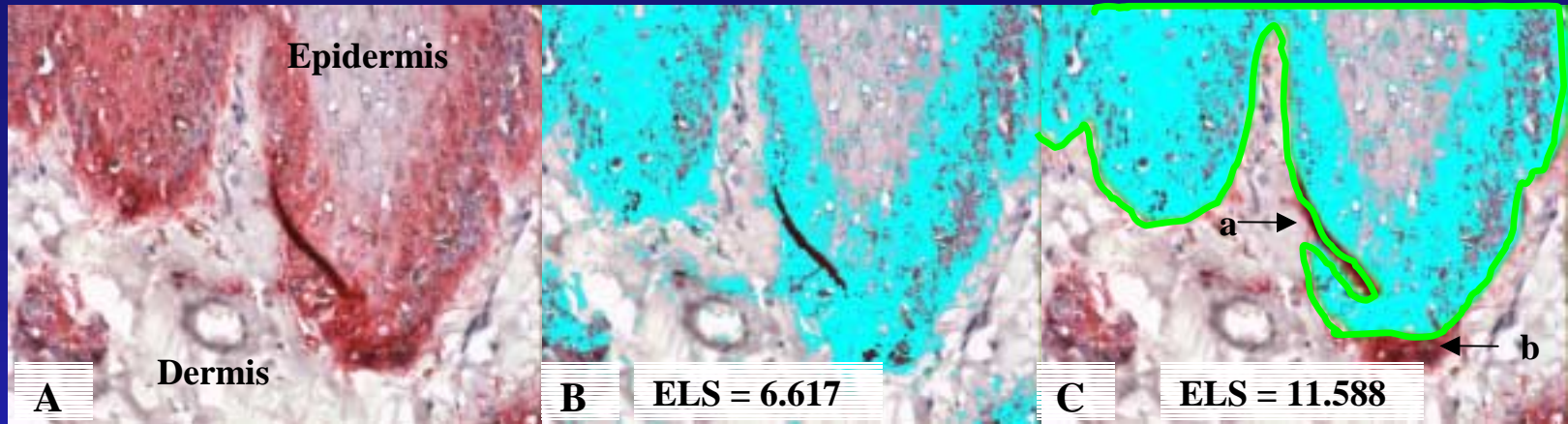


Figure 4. Quantification of metallothionein expression in *Stenella attenuata* skin using computerized image analysis, comparing non-outlined and outlined images.

The expression of MT was visualized as a red pigment in the immature keratinocytes in the epidermis (panel A). In panel B, the ELS was calculated from the signal in both the epidermis and the dermis. In panel C, the ELS was determined only from the MT signal within the outlined epidermal area-of interest. Outlining also allows artifacts of staining, which would alter the true expression level score, to be removed from the analysis. Such artifacts include folds in the skin, areas in which the skin has lifted up (a and b, panel C), tears, and edge effects. (Original magnifications: x400).

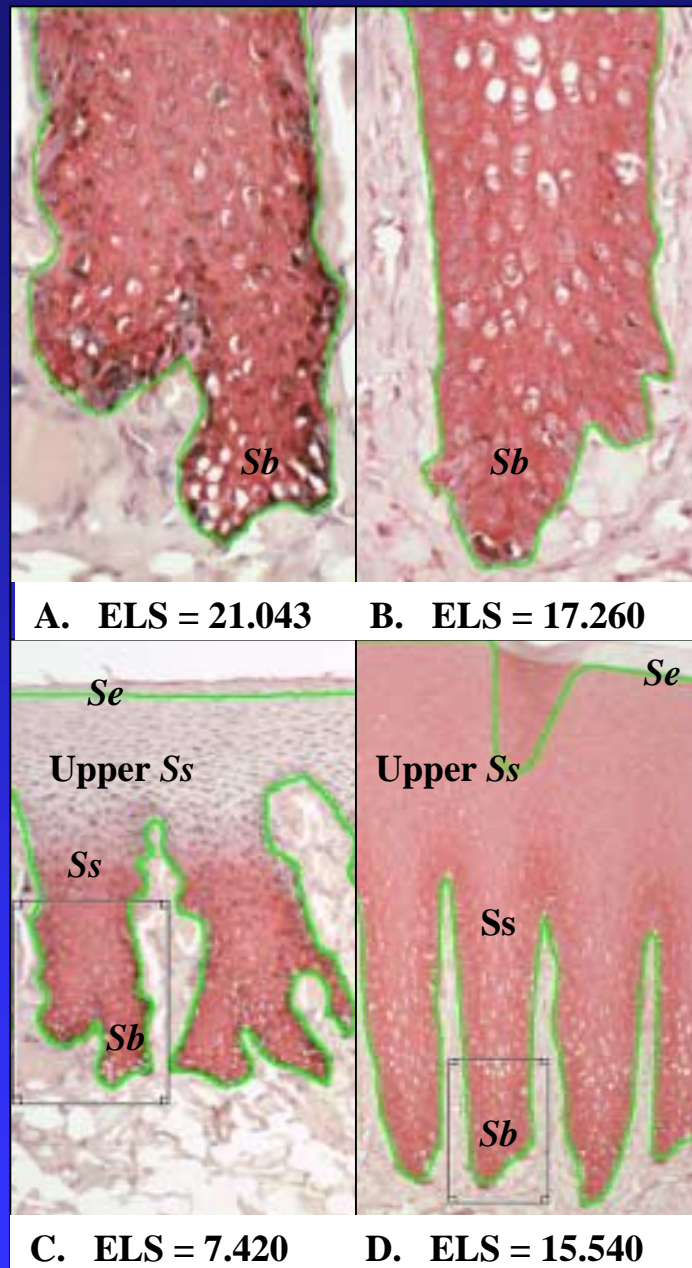


Figure 5. Comparison of CIA for two *Stenella attenuata* samples at different magnifications.

Many stress-response proteins can be expressed in all layers of the epidermis, from the *Stratum basale* (*Sb*), through the *Stratum spinosum* (*Ss*), to the *Stratum externum* (*Se*) (panel C). When two samples at 400x magnification were compared, similar ELS values were obtained for expression in the *Stratum basale* (panels A and B). The same samples analyzed at 100x show very different patterns of cellular accumulation of antigen, reflected by their much more representative ELS values (panels C and D). When staining is homogeneous and the cellular area of staining is the same on the sections being compared, a higher magnification can be used. When there is variability in the staining pattern, however, a lower magnification can discriminate staining differences better. (Samples stained with a combination of 40 SRP's).

Summary of Results

- Computer-assisted image analysis can be used to quantify immunohistochemical signal in cetacean skin.
- CIA results corresponded significantly with the subjective visual assessment ($0.001 < p < 0.01$), (**Figure 2**).
- Both the subjective visual assessment and the CIA methods accurately separated known healthy and diseased animals (**Figure 3**).
- ELS values obtained from the CIA method in two separate sessions were shown to be significantly concordant ($0.001 < p < 0.01$).
- This CIA method allows calculations to be performed within a selected area-of-interest, avoiding artifacts of staining which can confound results (**Figure 4**).
- CIA can be used at different magnifications, which is important for non-homogeneous stainings (**Figure 5**).

Applications

- At present, this technique of computer-assisted image analysis to quantify molecular stress has been applied to evaluate the impact of the tuna fishery on spotted dolphins in the Eastern Tropical Pacific, and the emergence of emaciation and low reproduction in the North Pacific gray whales.

ACKNOWLEDGEMENTS

Southwest Fisheries Science Center, NMFS, La Jolla, California

- Dr. Greg O'Corry-Crowe
- Dr. Bill Perrin
- Robin LeRoux
- Henry Orr
- Roy Allen

Media Cybernetics, Silver Spring, Maryland

- Will Casavan